

DENITROGENATION OF MODEL OIL USING IMIDAZALIUM-BASED IONIC LIQUID

MOHD HANAFI BIN SHAARI

**BACHELOR OF CHEMICAL ENGINEERING (GAS TECHNOLOGY)
UNIVERSITI MALAYSIA PAHANG**

©MOHD HANAFI BIN SHAARI (2013)

Thesis Access Form

No _____ Location: Universiti Malaysia Pahang

Author: Mohd Hanafi bin Shaari

Title: Denitrogenation of Model Oil Using Imidazalium-based Ionic Liquid

Status of access OPEN / RESTRICTED / CONFIDENTIAL

Moratorium period: _____ years, ending _____ / _____ 200_____

Conditions of access proved by (CAPITALS):

Supervisor (Signature).....

Faculty: Faculty of Chemical & Natural Resources Engineering (FKKSA)

Author's Declaration: *I agree the following conditions:*

OPEN access work shall be made available (in the University and externally) and reproduced as necessary at the discretion of the University Librarian or Head of Department. It may also be copied by the British Library in microfilm or other form for supply to requesting libraries or individuals, subject to an indication of intended use for non-publishing purposes in the following form, placed on the copy and on any covering document or label.

The statement itself shall apply to ALL copies:

This copy has been supplied on the understanding that it is copyright material and that no quotation from the thesis may be published without proper acknowledgement.

Restricted/confidential work: All access and any photocopying shall be strictly subject to written permission from the University Head of Department and any external sponsor, if any.

Author's signature.....**Date:**

users declaration: for signature during any Moratorium period (Not Open work):

I undertake to uphold the above conditions:

Date	Name (CAPITALS)	Signature	Address

DENITROGENATION OF MODEL OIL USING IMIDAZALIUM-BASED IONIC LIQUID

MOHD HANAFI BIN SHAARI

Thesis submitted in partial fulfilment of the requirements
for the award of the degree of
Bachelor of Chemical Engineering (Gas Technology)

Faculty of Chemical & Natural Resources Engineering

UNIVERSITI MALAYSIA PAHANG

JULY 2013

©MOHD HANAFI BIN SHAARI (2013)

SUPERVISOR'S DECLARATION

We hereby declare that we have checked this thesis and in our opinion, this thesis is adequate in terms of scope and quality for the award of the degree of Bachelor of Chemical Engineering (Gas Technology).

Signature :
Name of main supervisor : DR. SYAMSUL BAHARI ABDULLAH
Position :
Date :

STUDENT'S DECLARATION

I hereby declare that the work in this thesis is my own except for quotations and summaries which have been duly acknowledged. The thesis has not been accepted for any degree and is not concurrently submitted for award of other degree.

Signature :.....
Name : MOHD HANAFI BIN SHAARI
ID Number : KC09050
Date : JULY 2013

Dedication

This thesis is dedicated to my respective parents, supervisor and friends who have been my constant source of inspiration. They have given me the drive and discipline to tackle any task with enthusiasm and determination. Without their love and support this project would not have been made possible.

ACKNOWLEDGEMENT

I would like to thank the following people and organisations;

- My supervisor DR. Syamsul Bahari Abdullah for the support that really helped before, during and after doing the Undergraduate Research Project I and Undergraduate Research Project II.
- Lab Coordinator from chemical lab of Faculty of Chemical and Natural Resources that gave a good cooperation to me when I wanted to use and borrow the equipment in the lab.
- Faculty of Industry and Science Technology for giving me the permission to use Fourier Transform Infrared Spectroscopy (FTIR) to get the results to analyse the content of the Ionic Liquid, 3-Ethylimidazolium-Ethyl Phosphate.
- My friend Wan Mohd Syazwan bin Wan Abdullah for teaching me how to make the calculation for the sample for using in the High Performance Liquid Chromatography (HPLC) and to teach me on how to prepare the sample.

ABSTRACT

Removal of aromatic nitrogen contaminants from petroleum is important for many reasons. First, their combustion leads directly to the formation of nitrogen oxides (NO_x); emissions of NO_x, which contributes to acid rain, are under increasingly stringent control by environmental regulation. This thesis contains the study of Model Oil (iso-octane) by extraction with imidazolium-based ionic liquid. The objective of this experiment is to determine the nitrogen compound (Pyridine) removal efficiency of imidazolium-based ionic liquid (3-ethylimidazolium-ethyl phosphate) from iso-octane. In order to characterize the imidazolium based ionic liquids, the process of characterization of 3-ethylimidazolium-ethyl phosphate was conducted by using Fourier Transform Infra Red (FTIR) spectroscopy.

ABSTRAK

Penyingkiran nitrogen aromatik daripada petroleum adalah penting untuk pelbagai sebab. Pertama sekali, pembakarannya menuju terus kepada pembentukan nitrogen oksida; pembebasan nitrogen oksida menyumbang kepada hujan asid yang mana di kawal ketat oleh peraturan alam sekitar. Tesis ini mengandungi kajian mengenai model minyak (iso-oktana) melalui pengekstrakan dengan cecair ionik teras imidazalium. Objektif eksperimen ini adalah untuk menentukan kecekapan penyingkiran kompaun nitrogen (pyridine) daripada cecair ionik teras imidazalium (3-etilimidazalium-etil fosfat) daripada iso-oktana. Untuk mencirikan cecair ionik teras imidazalium, proses mencirikan 3-etilimidazalium-etil fosfat adalah di jalankan menggunakan Fourier Transform Infra Red (FTIR) spectroscopy.

TABLE OF CONTENTS

SUPERVISOR'S DECLARATION.....	IV
STUDENTS DECLARATION.....	V
<i>Dedication</i>	VI
ACKNOWLEDGEMENT.....	VII
ABSTRACT.....	VIII
ABSTRAK.....	IX
TABLE OF CONTENTS.....	1
LIST OF FIGURES.....	2
LIST OF TABLES.....	3
LIST OF ABBREVIATIONS.....	4
1. INTRODUCTION.....	5
1.1. Motivation and statement of problem.....	5
1.2. Objectives.....	6
1.3. Scope of study.....	6
1.4. Main contribution of this work.....	6
1.5. Organisation of this thesis.....	6-7
2. LITERATURE REVIEW.....	8
2.1 Overview.....	8
2.2 Introduction.....	8
2.3 Previous work on denitrogenation.....	9-10
2.4 Summary.....	11
3. MATERIALS AND METHODS.....	12
3.1 Overview.....	12
3.2 Introduction.....	12-14
3.3 Chemicals.....	14-16
3.4 Summary.....	17
4. DENITROGENATION OF MODEL OIL.....	18
4.1 Overview.....	18
4.2 Introduction.....	19
4.2.1 Results of FTIR spectroscopy analysis.....	19-20
4.3 Summary.....	21
5. CONCLUSION.....	22
5.1 Conclusion.....	22
5.2 Future work.....	22
REFERENCES.....	23-25
APPENDICES.....	26-30

LIST OF FIGURES

Figure 2-1: Extraction of the model oil containing indole as the N-compound with IL ([C2 mim] [Cl]) at 60°C.....	9
Figure 2-2: Extraction of the model oil containing pyridine as the N-compound with IL ([C2 mim] [Cl]) at 60°C.....	10
Figure 3-1: Calibration curve from HPLC results.....	14
Figure 4-1: FTIR analysis result.....	19
Figure 4-2: Structure of 3-ethylimidazolium-ethyl phosphate.....	20

LIST OF TABLES

Table 4-1: Analysis of FTIR spectroscopy result.....	20
--	----

LIST OF ABBREVIATIONS

C4	Butane
NO _x	Nitrogen oxides
HPLC	High Performance Liquid Chromatography
FTIR	Fourier Transform Infrared
MALDI-TOF	Matrix-assisted laser desorption-ionization time-of-flight mass spectrometry
GC	Gas Chromatography
IR	Infrared
C	Carbon
H	Hydrogen
N	Nitrogen
IL	Ionic liquid
S	Sulphur
wt-ratio	weight ratio
ppm	part pe million
M1V1	concentration 1×volume 1

1.0 INTRODUCTION

1.1 Motivation and statement of problem

Developments in the field of separation of pollutants, such as nitrogen and sulphur-containing compounds in the refinery industry have been experiencing a rapid growth due to more and stricter environmental protection regulations on the release of these pollutants from transportation fuels. Ionic liquids on the other hand are potential alternative solvents for the separation of aromatic hydrocarbons (benzene, toluene, ethyl benzene) from C4 to C10 aliphatic hydrocarbon mixtures. The prospect of ionic liquids as alternative solvents to substitute less desirable organic media is under intensive research. The low vapor pressure and ability to adjust the physical and chemical properties of ionic liquids by selection of the cations and anions, makes them attractive for a wide range of applications. For example, ionic liquid have already illustrated their potential as alternative solvents and reaction media. Furthermore, ionic liquids have been applied in solvent extraction processes and in gas separations (Ikenna Anugwom *et al.* 2011).

In fuels, nitrogen compounds act as precursors for nitrogen oxides (NO_x) which are environmental pollutants. Nitrogen are responsible for the formation of smog, sour gases, acid rain and NO_x emissions, these nitrogen compounds needs to be removed. Ionic liquid can be used to absorbed nitrogen compounds from fuels. Therefore, in this study is to investigate the capability of imidazolium-based ionic liquid in removing nitrogen compound from fuel.

1.2 Objectives

The following are the objectives of this research:

- To investigate the potential of imidazolium-based ionic liquids in removing nitrogen compound from fuel.

1.3 Scope of Study

The following are the scope of this research:

- i) To characterize the synthesized ionic liquids, 3-Ethylimidazolium-Ethyl Phosphate.
- ii) To study the removal performance of nitrogen compound from Iso-Octane.

1.4 Main contribution of this work

This experiment main purpose is to find the right imidazolium-based ionic liquid that will remove the some percentage of nitrogen compound from fuel. This experiment use 3-ethylimidazolium-ethyl phosphate, we will test with the fuel that contain nitrogen compound. The percentage of the decrease or increase amount of nitrogen in the mixture will be check by refer to HPLC result. From the graph we get from HPLC data results,we can know the outcome of this experiment.

1.5 Organisation of this thesis

The structure of the reminder of the thesis is outlined as follow:

Chapter 2 provides a description of review about ionic liquid characteristic and its usage. Also review about the function of FTIR and HPLC spectroscopy. Previous work for denitrogenation also state in this chapter. Use imidazolium-based ionic liquid too. The figure of the result also being show in the previous denitrogenation experiment.

Chapter 3 gives a review of materials and methods use for this experiment. Material use are by borrow from the chemical lab store. Methods involve are caharacterization of ionic liquid, preparation of standard sample, preparation of sample of model oil + pyridine + ionic liquid,

hplc spectroscopy and preparation of mobile phase. Chemicals use in this experiment are iso-octane, pyridine, 3-ethylimidazolium-ethyl phosphate and methanol.

Chapter 4 is review about Denitrogenation of model oil. This chapter focus on the results of this experiment. Consist of results of FTIR spectroscopy spectroscopy. Also contain the possible structure of 3-ethylimidazolium-ethyl phosphate and the table of wavelength and their respective molecular motion and functional group.

Chapter 5 give a review about conclusion of the experiment base on the result get in the chapter four. Consist of conclusion and future work recommendation step to make the experiment more succesful and accurate.

2.0 LITERATURE REVIEW

2.1 Overview

Ionic Liquid is a salt that are poorly coordinated, which cause it become liquid if temperature is below 100°C or even at room temperature. If at least an ionic liquid contain one ion that has delocalized charge and organic- can prevents formation of a stable crystal lattice to ionic liquid. ionic liquid , example 1-ethyl-3-methylimidazolium ethylsulphate-melting point <20 °C. Inorganic salt like sodium chloride have melting point at 801 °C. Instead of melting point a glass transition is observed especially if long aliphatic side chains are involved. Special characteristics of ionic liquid: thermal stability, low vapor pressure, electric conductivity, interesting solvent properties, liquid crystalline structures, high electroelasticity, high heat capacity and non flammability. Applications of ionic liquid: for lubricants and additive, use for lubricants and fuel additives, for electroelastic material, use for artificial muscles and robotics, for analytics, use for MALDI-TOF-matrices, GC-head-space-solvents and protein crystallization, for liquid crystals, use for displays, for heat storage, use for thermal, for solvents, use for Bio-catalysis, organic reaction and catalysis, Nano particles synthesis and polymerization, fluids for electrolytes, use for fuel cells, sensors, batteries, supercaps, metal finishing and coating, for separation, use for gas separations, extractive distillation, extraction and membranes. (ALDRICH, 2005)

2.2 Introduction

This paper present the result of the FTIR spectroscopy. From FTIR result, they are show the graph of different wavelength for different element. The element for each wavelength can be refer by table of IR Absorptions for Representative Functional Groups. This table contains info about group of alkanes, alkenes, alkynes, aromatics, alcohols, ethers, aldehydes, ketones, carboxylic acids, esters, acid chlorides, anhydrides, amines, amides, alkyl halides, nitriles, isocyanates, isothiocyanates, imines, nitro groups, mercaptans, sulfoxides, sulfones, sulfonates, phosphines and phosphine oxides. This table shows the molecular motion (etc: C-H stretch, =CH stretch) and wavenumber in unit of cm^{-1} .

2.3 Previous Work on Denitrogenation

Liquid-liquid extraction by using 1-ethyl-3-methylimidazolium chloride [C2 mim] [Cl] was found to be a very promising method for the removal of N- and S-compounds. This was evaluated by using a model oil (dodecane) with indole as a neutral nitrogen compound and pyridine as a basic nitrogen compound. An extraction capacity of up to 90 wt% was achieved for the model oil containing pyridine. The presence of aromatic compounds in the model oil decreases the extraction capacity of the ionic liquid. The presence of an aromatic compound in the sample affected the extraction capacity of the ionic liquid probably due to the fact that the ionic liquid exhibits a high affinity towards the aromatics. The effect of the ionic liquid-to-oil ratio was investigated in the extraction of N-compounds from the model oil. The extraction capacity improved as the ratio of IL-to-oil was increased, for indole Figure 1 and pyridine Figure 2, respectively. The extraction capacity also improved with an increase in the IL to oil ratio.

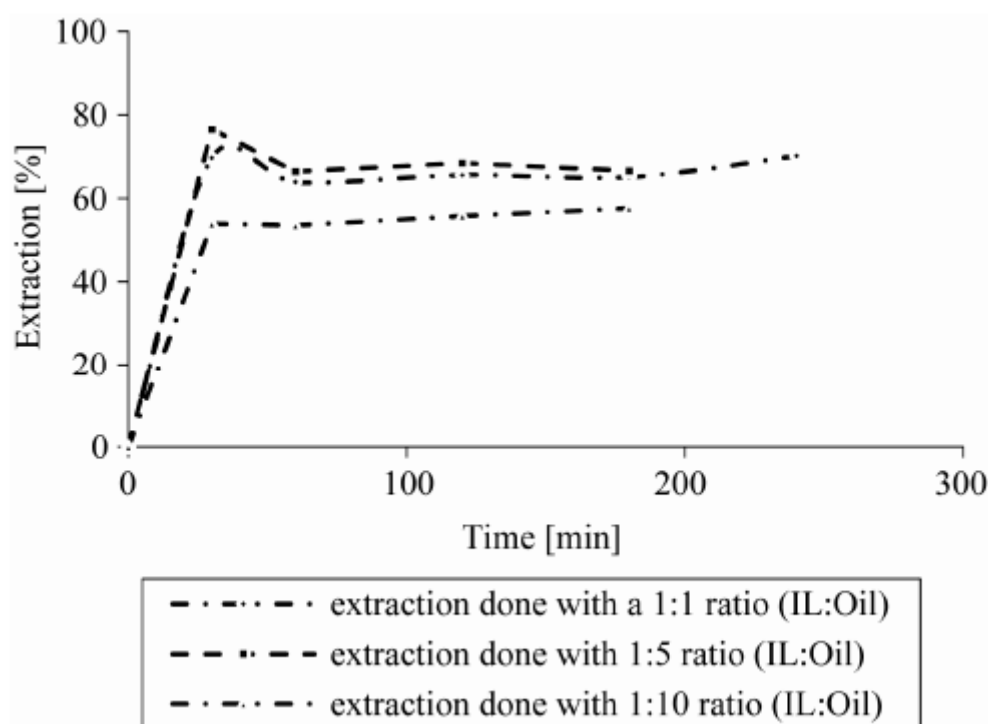


Figure 2-1: Extraction of the model oil containing indole as the N-compound with IL ([C2 mim] [Cl]) at 60°C.

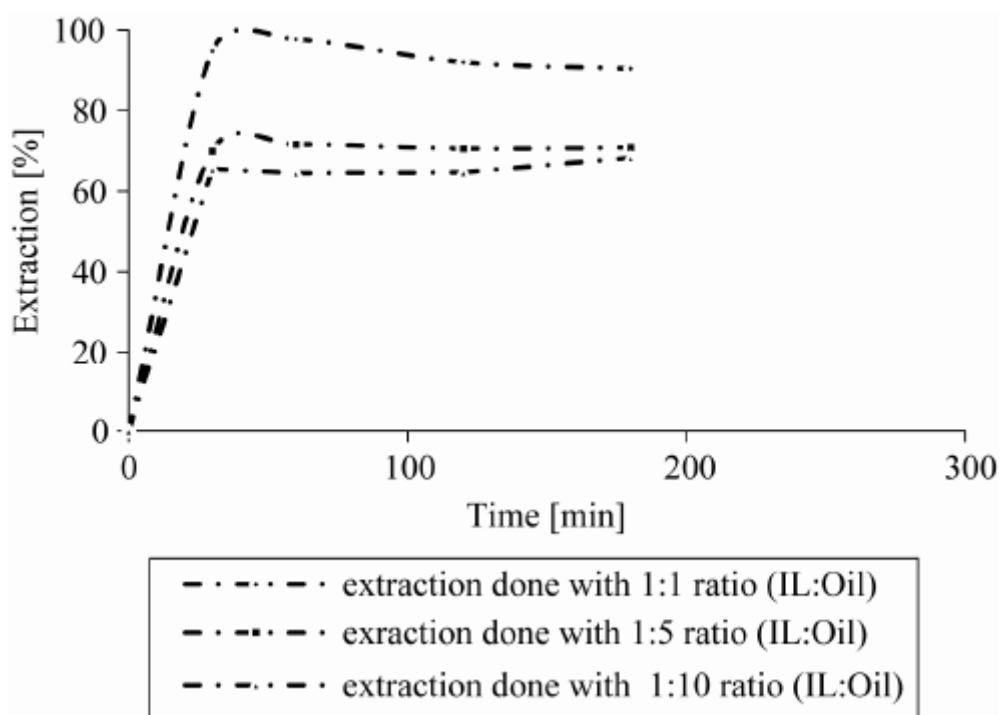


Figure 2-2: Extraction of the model oil containing pyridine as the N-compound with IL ([C2 mim] [Cl]) at 60°C.

The model oil in the current study was a mixture containing 3000 ppm of the nitrogen compound (indole or pyridine) in dodecane. Different ionic liquid to oil ratios were selected 1:1, 1:5, 1:10 (wt-ratio) and the extraction conditions were allowed to be the same in all cases: temperature was 60°C while 15 min mixing contact time between the oil and IL was used (Figures 1 and 2). The extraction of the model oil containing indole as the nitrogen compound was performed at 60°C and the IL-to oil ratio (wt) was 1:5. In equilibrium, 70% of the nitrogen compound was removed (Figure 1). In the case of a 1:1 IL-to-oil ratio 76% of indole was removed in equilibrium at 60°C. This is an indication that a higher IL-to-oil ratio does not improve the extraction effectiveness prominently. Moreover, thinking about potential large-scale process, 1:5 IL to oil ratio is already rather high.

2.4 Summary

Ionic Liquid is a salt that are poorly coordinated, which cause it become liquid if temperature is below 100°C or even at room temperature. Applications of ionic liquid: for lubricants and additive, for electroelastic material, for analytics, for liquid crystals, for solvents, fluids for electrolytes and for gas separations. FTIR shows the elemental group occur in the ionic liquid use and HPLC shows whether the pyridine (nitrogen) is reduce or not from the graph that have been constructed. From the theory it looks like extraction capacity improved as the ratio of IL-to-oil was increased but from the previous experiment result a higher IL-to-oil ratio does not improve the extraction effectiveness prominently.

3.0 MATERIALS AND METHODS

3.1 Overview

The 3-ethylimidazolium-ethyl phosphate is analyzed by the FTIR spectroscopy to know the elemental groups that occur in that ionic liquid.

3.2 Introduction

3.2.1 Characterization of Ionic Liquid

Some sample of the 3-ethylimidazolium was put on the sample plate in the FTIR equipment. Close the cover of the place of the sample plate. From the computer click start button to begin the simulation. As the graph of wavelength occurs, the data was saved in the compact disk for further analysis later. The analysis of the wavelength data can be done by referring to the table of IR Absorptions for Representative Functional Groups. The functional group and molecular motion were found based on the table of IR Absorptions for Representative Functional Groups.

3.2.2 Preparation of standard sample

The pyridine needed to be add to the 100mL of model oil to get 1500ppm was calculated by using formula:

$$1500ppm = \frac{x \text{ mg pyridine}}{0.1\text{mL model oil}}$$

Conversion: density = mass/volume

Volume pyridine = mass pyridine/density pyridine

The 100mL mixture was divided into two container for each container contain 50mL. One of the 50mL mixture was used to make concentration of 1000ppm, 800ppm, 500ppm and 250ppm while the other one 50mL mixture is save. By using formula $M_1V_1 = M_2V_2$, the V_2 needed was got to produce 1000ppm, 800ppm, 500ppm and 250ppm concentration of model oil mixture by assume V_1 for some value. V_2 is the total volume needed to produce wanted concentration model oil mixture by adding it with 1500ppm model oil mixture for amount of V_1 we assumed earlier. $V_2 - V_1$ is equal to volume needed more from 1500ppm model oil mixture to produce certain concentration. After all the standard sample are prepared, five sample from the five different concentration were made. All preparation of sample into vial bottles must be made by using filter to avoid any the failure of the HPLC equipment.

3.2.3 Preparation of sample of model oil + pyridine + ionic liquid

From sample prepared in the schott bottle 50mL, take volume needed to produce 1(ionic liquid):1(model oil mixture) ratio sample. Assume for 1 : 0.5mL. Prepare for one sample for ratio 1:1 for each concentration. The step 1 to 2 are repeated for the ratio of 1:4. After all the sample are prepared, the sample are place in the vials bottle by using syringe and filter.

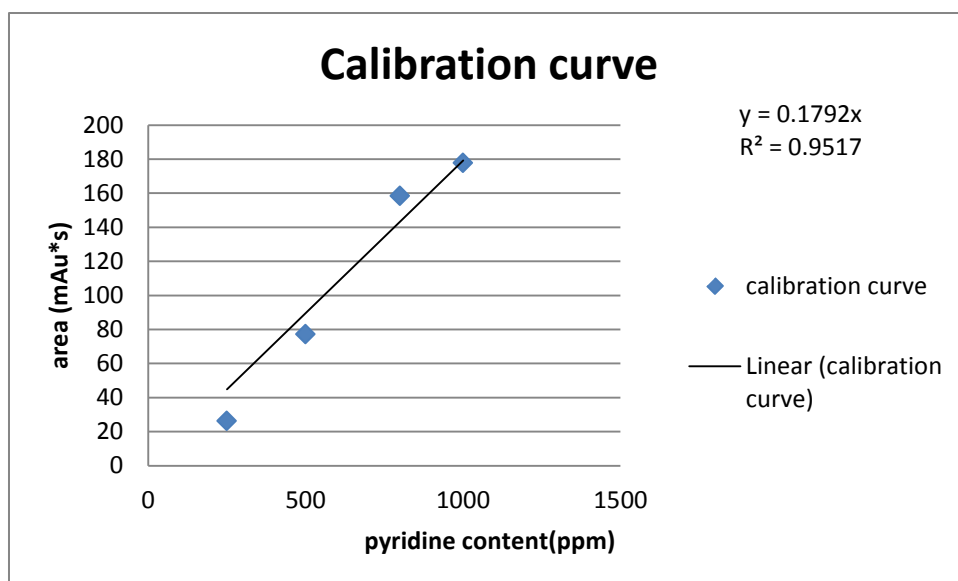


Figure 3-1: Calibration curve from HPLC results.

3.2.4 HPLC Spectroscopy

Ten sample were placed in the place of the vials bottle. They were placed according to the name starting from four sample of standard sample, four sample of ratio 1:1 and four sample of ratio 1:4. The HPLC spectroscopy equipment was run by lab coordinator. The results of the analysis of the ten sample were got after 24 hours of analysis.

3.2.5 Preparation of mobile phase

Mobile phase for this experiment was methanol. In the ratio of 1:9, 90% percent of methanol was used to produce the solution of mobile phase while other 10% ratio for ultra pure water. For produce 1 liter mobile phase, 100mL of ultra pure water was added into 900mL of methanol. The mobile phase solution was filtered by using membrane filter. Then the mobile phase solution was degas with sonicator.

3.3 Chemicals

3.3.1 Model oil

The model oil use for this experiment was iso-octane. It is flammable in liquid and vapor state. It can causes causes eye, skin, and respiratory tract irritation. Its vapors may cause drowsiness and dizziness. Can cause aspiration hazard if swallowed as it can enter lungs and cause damage. It is very toxic to aquatic organisms as it may cause long-term adverse effects in the aquatic environment. For handling, only use it under a chemical fume hood together

with explosion-proof equipment. The user should wear personal protective equipment. It should keep away from open flames, hot surfaces and sources of ignition. For storage of iso-octane, it must keep away from open flames, hot surfaces and sources of ignition. The containers of iso-octane must tightly closed in a dry, cool and well-ventilated place.

3.3.2 Pyridine

Pyridine was used as source of nitrogen in this experiment. It is flammable in liquid and vapor state. Its vapor may cause flash fire. If inhale it can cause irritation, headache, drowsiness, dizziness and loss of coordination for short term exposure and for long term exposure can cause nausea, vomiting, diarrhea, stomach pain, loss of appetite, dizziness, sleep disturbances, emotional disturbances, loss of coordination and nerve damage. If the effects still occur, the victim should move to uncontaminated area. If breathing is difficult, oxygen should be administered by qualified personnel. The victims should get immediate medical attention.

3.3.3 3-Ethylimidazolium-ethyl phosphate

Ionic liquid use for this experiment was 3-ethylimidazolium-ethyl phosphate. It is corrosive to skin, eyes, and respiratory system. Its liquid or spray mist may produce tissue damage, particularly in mucous membranes of the eyes, mouth and respiratory tract. If contact to the skin may produce burns. If contact to the eyes can result in corneal damage or blindness. Inhalation of the spray mist may produce severe irritation of respiratory tract, characterized by coughing, choking, or shortness of breath. It is a corrosive materials that may cause serious injury if ingested. If contact to the eyes occur, remove any contact lenses. Flush eyes with running water for a minimum of 15 minutes, occasionally lifting the upper and lower eyelids. If the chemical gets spilled on a clothed portion of the body, remove the contaminated clothes as quickly as possible and protecting your own hands and body. Place the victim under a deluge shower. If the chemical touches the victim's exposed skin, such as the hands: Gently and thoroughly wash the contaminated skin with running water and non-abrasive soap.

3.3.4 Methanol

A mobile phase used in the HPLC equipment is methanol. It should not use in a confined area without proper ventilation. If wear the contact lenses it may cause further damage in case of splash into eye. It should avoid use near heat, flames, sparks, and other sources of ignition. It is extremely flammable in liquid and vapour state. It is fatal if swallowed. It may damage fertility or the unborn child (fetotoxic and teratogenic effects), cause damage to eyes and central nervous system if ingested or inhaled. If contact with eye, remove contact lenses if worn. In case of contact, immediately flush eyes with plenty of clean running water for at least 15 minutes, lifting the upper and lower eyelids occasionally. If contact with skin, remove contaminated clothing. In a shower, wash affected areas with soap and water for at least 15 minutes. Seek medical attention if irritation occurs or persists. Wash clothing before reuse.

3.4 Summary

In order to know the percent of reduce nitrogen compound in the model oil it is needed to prepared standard sample and two different ratio of ionic liquid + model oil mixture. From the standard sample it will become the references to know wether the nitrogen removal was succeeded or not. The step involve for this experiment were characterization of ionic liquid by FTIR spectroscopy, preparation for different concentration of ionic liquid, preparation of some sample of different ratio and the analysis of the samples wtih HPLC spectroscopy. This experiment was done by using the model oil,pyridine and ionic liquid.

4.0 Denitrogenation of Model Oil

4.1 Overview

This chapter review about the results of this thesis. The result get from the FTIR spectroscopy analysis. Result of FTIR analysis show the graph of different wavelength for different element. Its analyze the element wavelength in range of 700cm^{-1} to 4000cm^{-1} . The lowest wavelength is 748.11cm^{-1} which refer to aromatics, C-H bend (ortho) and the highest wavelength is 3391.50cm^{-1} which refer to alcohols, O-H stretch.

4.2 Introduction

4.2.1 Results of FTIR spectroscopy analysis

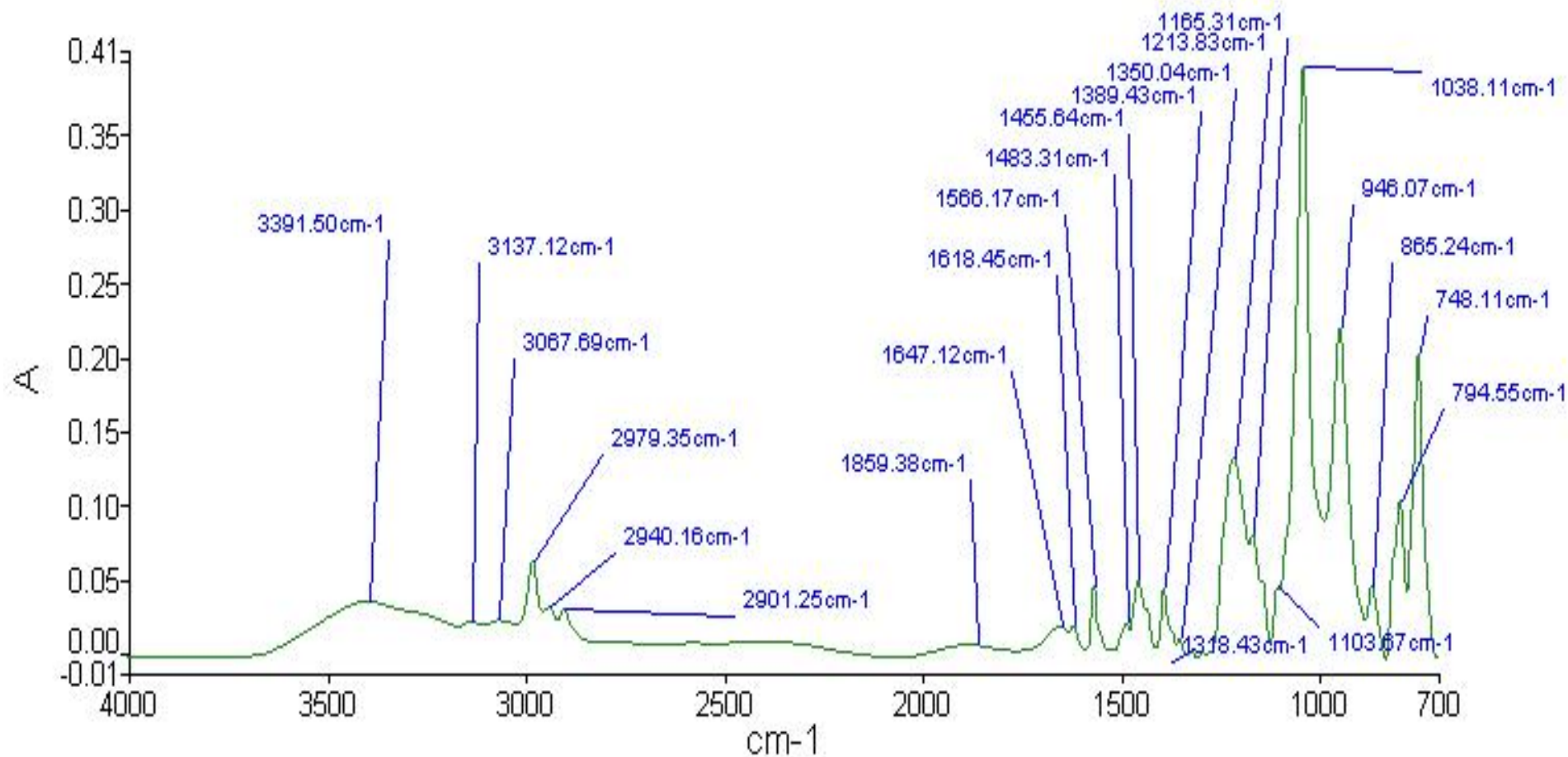


Figure 4-1: FTIR analysis result

Possible structure of 3-ethylimidazolium-ethyl phosphate

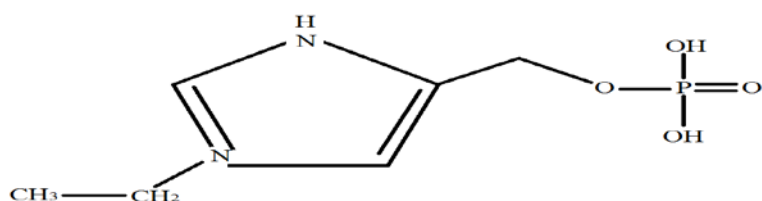


Figure 4-2: Structure of 3-ethylimidazolium-ethyl phosphate

Table 4-1: Analysis of FTIR spectroscopy result

Peak number	Wavenumber (cm ⁻¹)	Functional group	Molecular motion
1	748.11	aromatics	C-H bends (ortho)
2	794.55	amines	N-H bends
4	946.07	anhydrides	C-O stretch
6	1103.67	ketones	C-C stretch
7	1165.31	phosphine oxides	P=O
9	1318.43	amines	C-N stretch (aryl)
23	3391.5	alcohols	O-H

From the Table 4-1 analysis, it is known that from the possible structure of 3-ethylimidazolium-ethyl phosphate, this ionic liquid contains at wavelength 748.11cm⁻¹ its refer to C-H bend (ortho) from aromatics group, at wavelength 794.55cm⁻¹ its refer to N-H bend from amines group, at wavelength 946.07cm⁻¹ its refer to C-O stretch from anhydrides group, at wavelength 1103.67cm⁻¹ its refer to C-C stretch from ketones group, at wavelength 1165.31cm⁻¹ its refer to P=O from phosphine oxides group, at wavelength 1318.43cm⁻¹ its refer to C-N stretch (aryl) from amines group, and at wavelength 3391.5 its refer to O-H stretch from alcohols group.

4.3 Summary

Result of FTIR analysis show the graph of different wavelength for different element. Its analyze the element wavelength in range of 700cm^{-1} to 4000cm^{-1} . It is know that from possible structure of 3-ethylimidazolium-ethyl phosphate, this ionic liquid contain at wavelength 748.11cm^{-1} its refer to C-H bend (ortho) from aromatics group, at wavelength 794.55cm^{-1} its refer to N-H bend from amines group, at wavelength 946.07cm^{-1} its refer to C-O stretch from anhydrides group, at wavelength 1103.67cm^{-1} its refer to C-C stretch from ketones group, at wavelength 1165.31cm^{-1} its refer to P=O from phosphine oxides group, at wavelength 1318.43cm^{-1} its refer to C-N stretch (aryl) from amines group, and at wavelength 3391.5 its refer to O-H stretch from alcohols group.

5.0 CONCLUSION

5.1 Conclusion

3-ethylimidazolium-ethyl phosphate, this ionic liquid contain at wavelength 748.11cm^{-1} its refer to C-H bend (ortho) from aromatics group, at wavelength 794.55cm^{-1} its refer to N-H bend from amines group, at wavelength 946.07cm^{-1} its refer to C-O stretch from anhydrides group, at wavelength 1103.67cm^{-1} its refer to C-C stretch from ketones group, at wavelength 1165.31cm^{-1} its refer to P=O from phosphine oxides group, at wavelength 1318.43cm^{-1} its refer to C-N stretch (aryl) from amines group, and at wavelength 3391.5 its refer to O-H stretch from alcohols group.

5.2 Future Work

The volume of Ionic Liquid prepared should excess in volume so that experiment can be repeat more than two times if theres an error in the first experiment. The optimal absorption and desorption temperature are 298 K and 363 K, respectively. Make sure the temperature while using magnetic stirrer is constant at 323 K. Make sure the tube flow is good to avoid the ionic liquid absorb only little NO_2 in fuel. Ionic liquid must be dry for 2 days at 323 K to get the good ionic liquid solution mixture.

REFERENCES

1. Linzhi Zhai, Qin Zhong, Chuan He, Juan Wang- Hydroxyl ammonium ionic liquids synthesized by water-bath microwave: Synthesis and desulfurization, *Journal of hazardous material*. 177(2010) 807-813
2. P. Wasserscheid, W. Keim, *Angew-Review on Ionic Liquid, Chem. Int. Ed. Engl.* 39(2000) 1
3. Shuguang Zhang, Qinglin Zhang and Z. Conrad Zhang-Extractive Desulfurization and Denitrogenation of Fuels Using Ionic Liquids, *Ind. Eng. Chem. Res.* 43, (2004) 614-622
4. Aldrich Company USA-Enabling Technologies Ionic Liquid, *Chemical Synthesis*. Volume 5 No. 6 (2005) 1-24
5. I Wen Sun et. al-Cyclic Ammonium-Based Ionic Liquids as Potential Electrolytes for Dye-Sensitized Solar Cells, *Int. J. Electrochem. Sci.*, 7 (2012) 1-17
6. University of Wollongong (2003)-Ionic liquids dramatically improve performance of artificial muscles. Retrieved on 30 November 2012 from <http://media.uow.edu.au/releases/2002/artificial.html>
7. NDT Education Resource Center (2001-2012)-Electrical Conductivity and Resistivity. Retrieved on 30 November 2012 from http://www.ndt-ed.org/EducationResources/CommunityCollege/Materials/Physical_Chemical/Electrical.htm
8. Loughborough University- The Electrodeposition of Zinc-Manganese Alloy Coatings from Ionic Liquid Electrolytes 1-30
9. K. Gable (2000)-FTIR Spectroscopy. Retrieved on 30 November 2012 from <http://www.chemistry.oregonstate.edu/courses/ch361-464/ch362/irinstrs.htm>

10. Iolitec USA (2011)-Fuel Cells. Retrieved on 30 November 2012 from <http://www.iolitec-usa.com/Ionic-Liquids/fuel-cells.html>
11. George Mason University (1998)- Gas Chromatography-Mass Spectroscopy Background. Retrieved on 30 November 2012 from <http://www.gmu.edu/depts/SRIF/tutorial/gcd/gc-ms2.htm>
12. CMMAP (2000)-Heat Capacity. Retrieved on 30 November 2012 from <http://www.cmmap.org/learn/climate/heatCapacity.html>
13. Net Industries (2012)- Heat Capacity - Heat Capacity And Calorimetry, Heat Capacity And The Law Of Conservation Of Energy - Significance of the high heat capacity of water. Retrieved on 30 November 2012 from <http://science.jrank.org/pages/3265/Heat-Capacity.html>
14. Encyclopedia Britannica (2012)-Liquid Crystals. Retrieved on 30 November 2012 from <http://www.britannica.com/EBchecked/topic/343083/liquid-crystal>
15. Zeyun Wang et. al (2012)- Conductive Lubricating Grease Synthesized Using the Ionic Liquid, Tribol Lett (2012) 46:33–42 1-10
16. Gabriel Gache (2008)-Nitrogen Compound Pollution Can Accelerate Global Warming. Retrieved on 30 November 2012 from <http://news.softpedia.com/news/Nitrogen-Compounds-Pollution-Can-Accelerate-Global-Warming-85816.shtml>
17. Bruker Corporation (2012)-Biomolecular NMR. Retrieved on 30 November 2012 from <http://www.bruker.com/products/mr/nmr/dnp-nmr/overview.html>
18. Hampton Research (2003-2012)- Ionic Liquid Screen. Retrieved on 30 November 2012 from https://hamptonresearch.com/product_detail.aspx?cid=1&sid=198&pid=609

19. Ionic Liquid Technologies (2011)-Heat Transfer Media/Thermal Fluids. Retrieved on 30 November from <http://www.iolitec.de/en/Heat-Storage-Heat-Transport/heat-transfer-media-thermal-fluids.html>
20. Katherin Stillwater (1999-2012)-What can FTIR Spectroscopy Be Used For? Retrieved on 30 November 2012 from http://www.ehow.com/about_5659926_can-ftir-spectroscopy-used-for_.html
21. DR. Deepak Bhanot (2011-2012)-High Performance Liquid Chromatography: Module 3. Retrieved on 22 May 2013 from <http://lab-training.com/landing/free-hplc-training-programme-5/>
22. Ikenna Anugwom et. al (2011)-Ionic Liquid Assisted Extraction of Nitrogen and Sulphur-Containing Air Pollutants from Model Oil and Regeneration of the Spent Ionic Liquid, *Journal of Environmental Protection*, 2011, 2, 796-802.
23. Chemical Book (2008)-Ethyl Phosphate. Retrieved on 1 June 2013 from http://www.chemicalbook.com/ProductChemicalPropertiesCB0740449_EN.htm
24. Wikipedia (2013)-Ethyl Group. Retrieved on 1 June 2013 from http://en.wikipedia.org/wiki/Ethyl_group
25. Karolina Kedra Krolik et. al (2011)-Extraction of Thiophene or Pyridine from n-Heptane Using Ionic Liquids, Gasoline and Diesel Desulfurization, *Ind. Eng. Chem. Res.* 2011, 50, 2296–2306.
26. Thermo Nicolet Corporation (2001)-Introduction to Fourier Transform Infrared Spectrometry, 1-8.
27. Wikipedia (2013)-Imidazole. Retrieved on 1 June 2013 from <http://en.wikipedia.org/wiki/Imidazole>
28. Qilong Ren et. al (2011)- Ionic Liquid-Mediated Liquid-Liquid Extraction, 1-19.

APPENDICES

$$1. \quad 1500\text{ppm} = \frac{x \text{ mg pyridine}}{0.1 \text{ L Model oil}}$$

$$= 150 \text{ mg pyridine}$$

Density pyridine = mass pyridine/volume pyridine

$$150 \text{ mg pyridine} = 0.8046 \text{ mL pyridine}$$

$$2. \quad 1000\text{ppm}$$

$$M1V1 = M2V2$$

$$1500(0.020) = 1000V2$$

$$V2 = 0.03\text{L}$$

$$\text{Model oil} = 30\text{mL} - 20\text{mL} = 10\text{mL}$$

1500ppm: take 20mL

$$3. \quad 800\text{ppm}$$

$$M1V1 = M2V2$$

$$1500\text{ppm}(0.010) = 800V2$$

$$V2 = 0.019 = 0.020$$

$$\text{Model oil} = 20\text{mL} - 10\text{mL} = 10\text{mL}$$

1500ppm: take 10mL

$$4. \quad 500\text{ppm}$$

$$1500(0.008) = 500V2$$

$$V2 = 0.024$$

$$\text{Model oil} = 24\text{mL} - 8\text{mL} = 16\text{mL}$$

1500ppm: take 8mL

$$5. \quad 250\text{ppm}$$

$$1500\text{ppm}(0.004) = 250V2$$

$$V_2 = 0.024$$

$$\text{Model oil} = 24\text{mL} - 4\text{mL} = 20\text{mL}$$


1500ppm: take 4mL

$$\frac{x \text{ nitrogen}}{1500\text{ppm}} = \frac{14.0067 \text{ g/mol}}{79.0999 \text{ g/mol}}$$

$$x = 265.61 \text{ nitrogen}$$

PYRIDINE(PPM)	NITROGEN(PPM)
1500	265.61
1000	177.08
800	141.66
500	88.54
250	44.27

Results of FTIR

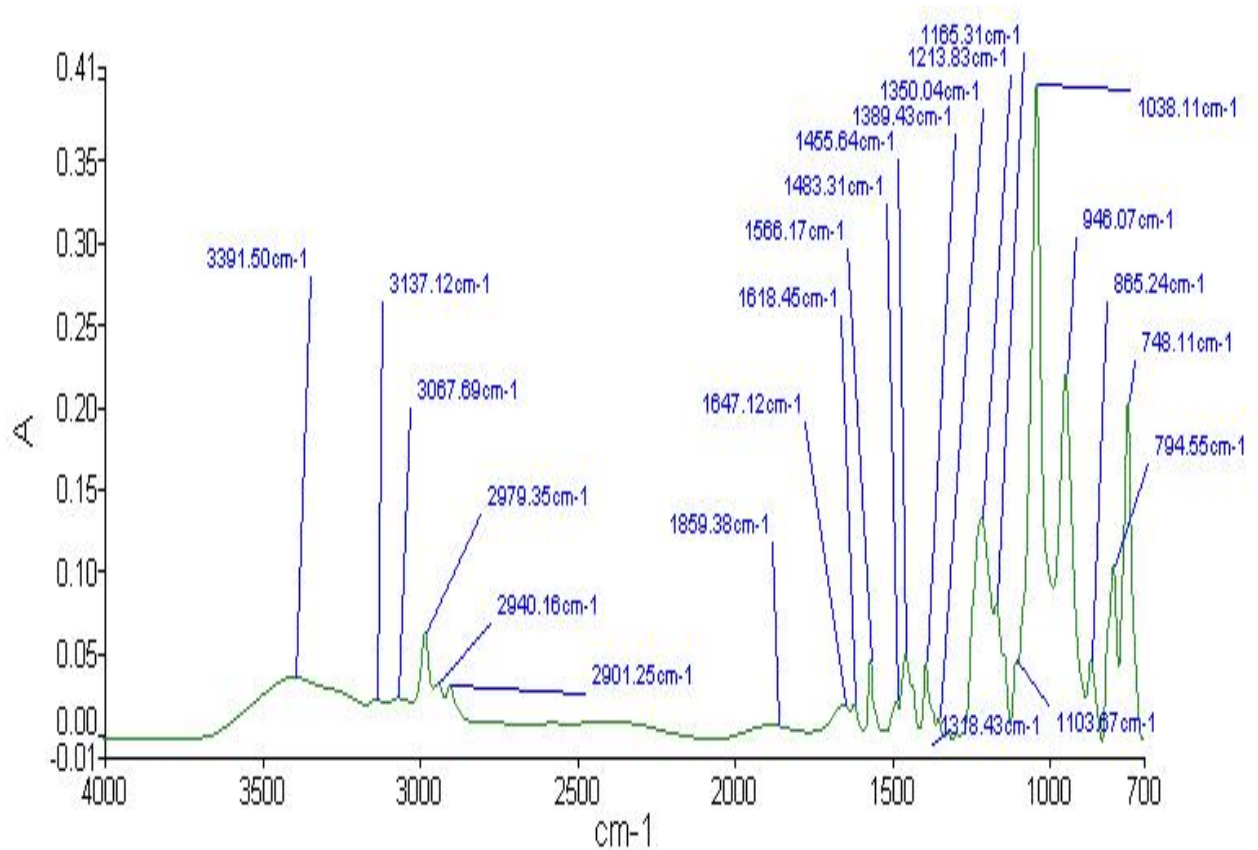

Infrared Spectroscopy

IR Absorptions for Representative Functional Groups

Functional Group	Molecular Motion	Wavenumber (cm ⁻¹)
alkanes	C-H stretch	2950-2800
	CH ₂ bend	~1465
	CH ₃ bend	~1375
	CH ₂ bend (4 or more)	~720
alkenes	=CH stretch	3100-3010
	C=C stretch (isolated)	1690-1630
	C=C stretch (conjugated)	1640-1610
	C-H in-plane bend	1430-1290
	C-H bend (monosubstituted)	~990 & ~910
	C-H bend (disubstituted - E)	~970
	C-H bend (disubstituted - 1,1)	~890
	C-H bend (disubstituted - Z)	~700
	C-H bend (trisubstituted)	~815
alkynes	acetylenic C-H stretch	~3300
	C,C triple bond stretch	~2150
	acetylenic C-H bend	650-600
aromatics	C-H stretch	3020-3000
	C=C stretch	~1600 & ~1475
	C-H bend (mono)	770-730 & 715-685
	C-H bend (ortho)	770-735
	C-H bend (meta)	~880 & ~780 & ~690
	C-H bend (para)	850-800
alcohols	O-H stretch	~3650 or 3400-3300
	C-O stretch	1260-1000
ethers	C-O-C stretch (dialkyl)	1300-1000

	C-O-C stretch (diaryl)	~1250 & ~1120
aldehydes	C-H aldehyde stretch	~2850 & ~2750
	C=O stretch	~1725
ketones	C=O stretch	~1715
	C-C stretch	1300-1100
carboxylic acids	O-H stretch	3400-2400
	C=O stretch	1730-1700
	C-O stretch	1320-1210
	O-H bend	1440-1400
esters	C=O stretch	1750-1735
	C-C(O)-C stretch (acetates)	1260-1230
	C-C(O)-C stretch (all others)	1210-1160
acid chlorides	C=O stretch	1810-1775
	C-Cl stretch	730-550
anhydrides	C=O stretch	1830-1800&1775-1740
	C-O stretch	1300-900
amines	N-H stretch (1 per N-H bond)	3500-3300
	N-H bend	1640-1500
	C-N Stretch (alkyl)	1200-1025
	C-N Stretch (aryl)	1360-1250
	N-H bend (oop)	~800
amides	N-H stretch	3500-3180
	C=O stretch	1680-1630
	N-H bend	1640-1550
	N-H bend (1°)	1570-1515
alkyl halides	C-F stretch	1400-1000
	C-Cl stretch	785-540
	C-Br stretch	650-510
	C-I stretch	600-485
nitriles	C,N triple bond stretch	~2250
isocyanates	-N=C=O stretch	~2270

isothiocyanates	-N=C=S stretch	~2125
imines	R ₂ C=N-R stretch	1690-1640
nitro groups	-NO ₂ (aliphatic)	1600-1530&1390-1300
	-NO ₂ (aromatic)	1550-1490&1355-1315
mercaptans	S-H stretch	~2550
sulfoxides	S=O stretch	~1050
sulfones	S=O stretch	~1300 & ~1150
sulfonates	S=O stretch	~1350 & ~11750
	S-O stretch	1000-750
phosphines	P-H stretch	2320-2270
	PH bend	1090-810
phosphine oxides	P=O	1210-1140